

MASTERflex[®]

Reimagine
Fluid Handling
Operate, collaborate and
connect with ease

Achrom



Masterflex[®] L/S[®]
MasterSense[™] Drives

Reimagine an easier way to pump with the **MasterSense drive**. Intuitive programming and forward-looking features enable you to work smarter and pump easier than ever before. And, with the same great features as all Masterflex peristaltic pump drives including high precision, repeatable flow control, and reliable performance, you can have it all in one small footprint—even more bench space.

Intuitive touch-screen interface

Easy programming functionality

Innovative forward-looking features

Multiple connectivity options

Versatile system components



FPO




Intuitive Touch-Screen Interface

- **Quick and easy** to set up and program
- 5" capacitive TFT/RGB touch screen with a clear, wide viewing angle
- Adjust parameters with a simple touch—**even with gloves on**

Easy & Accurate

- **Save and recall** up to 50 commonly used protocols, each with a unique name
- Operational modes include continuous run, volume dispense, and time dispense
- Batch control and interval/off time in volume and time modes
- Antidrip and speed ramping available in volume mode—ensure accuracy and **minimize sample loss**
- Know when you hit your mark—graphical indicator lets you know when loaded tubing is **optimally calibrated**


Continuous Run

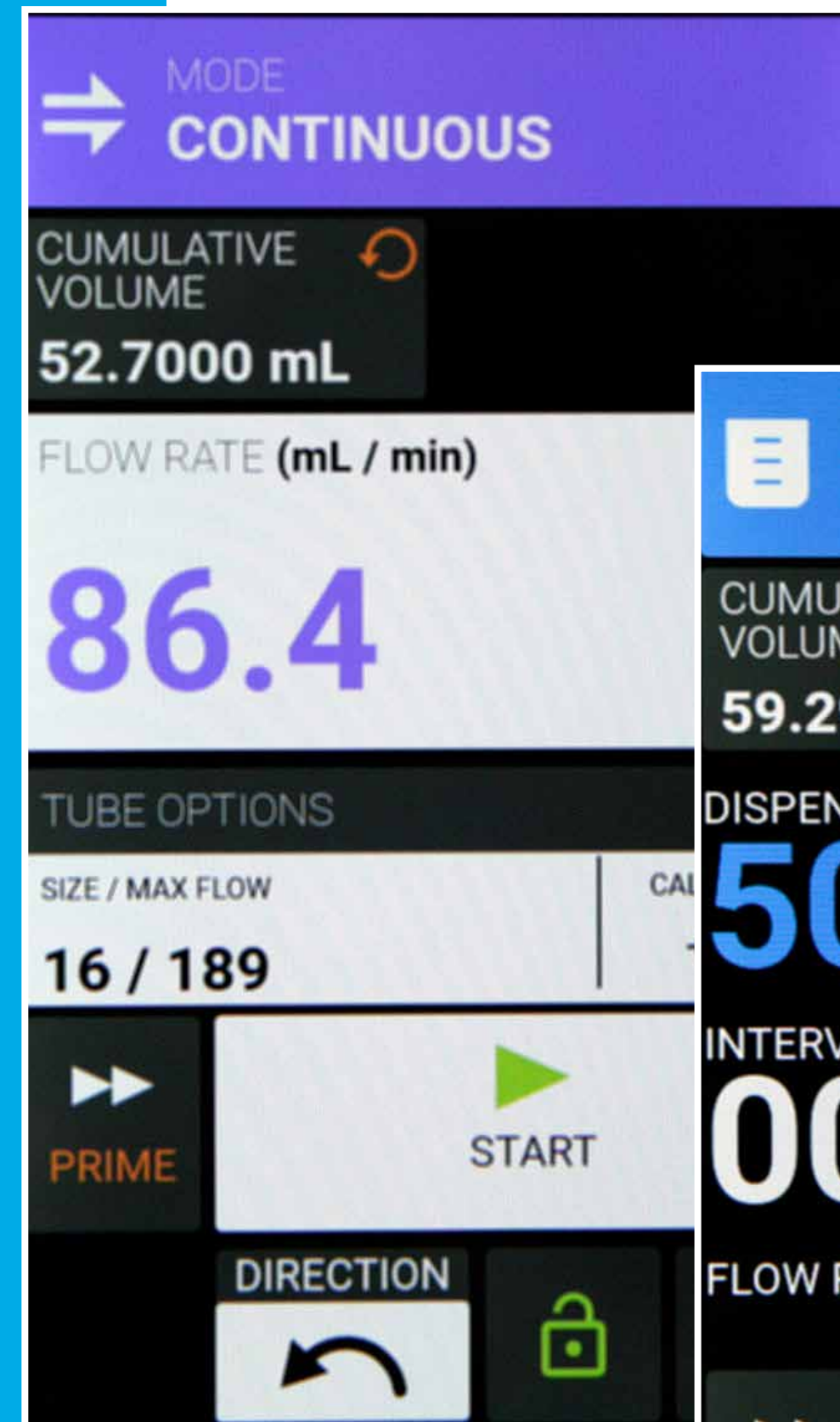


Time Dispense




Volume Dispense

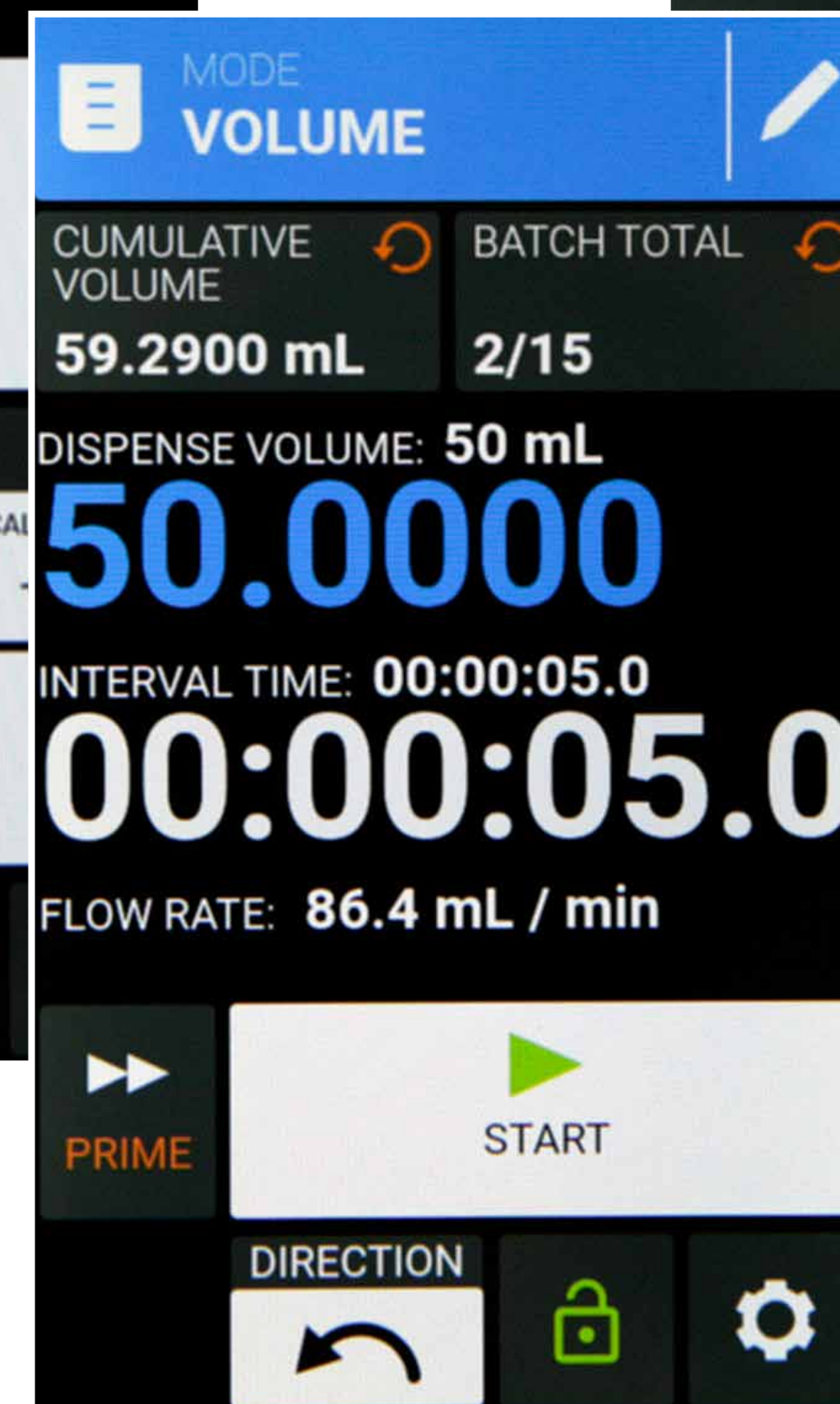

Calibration

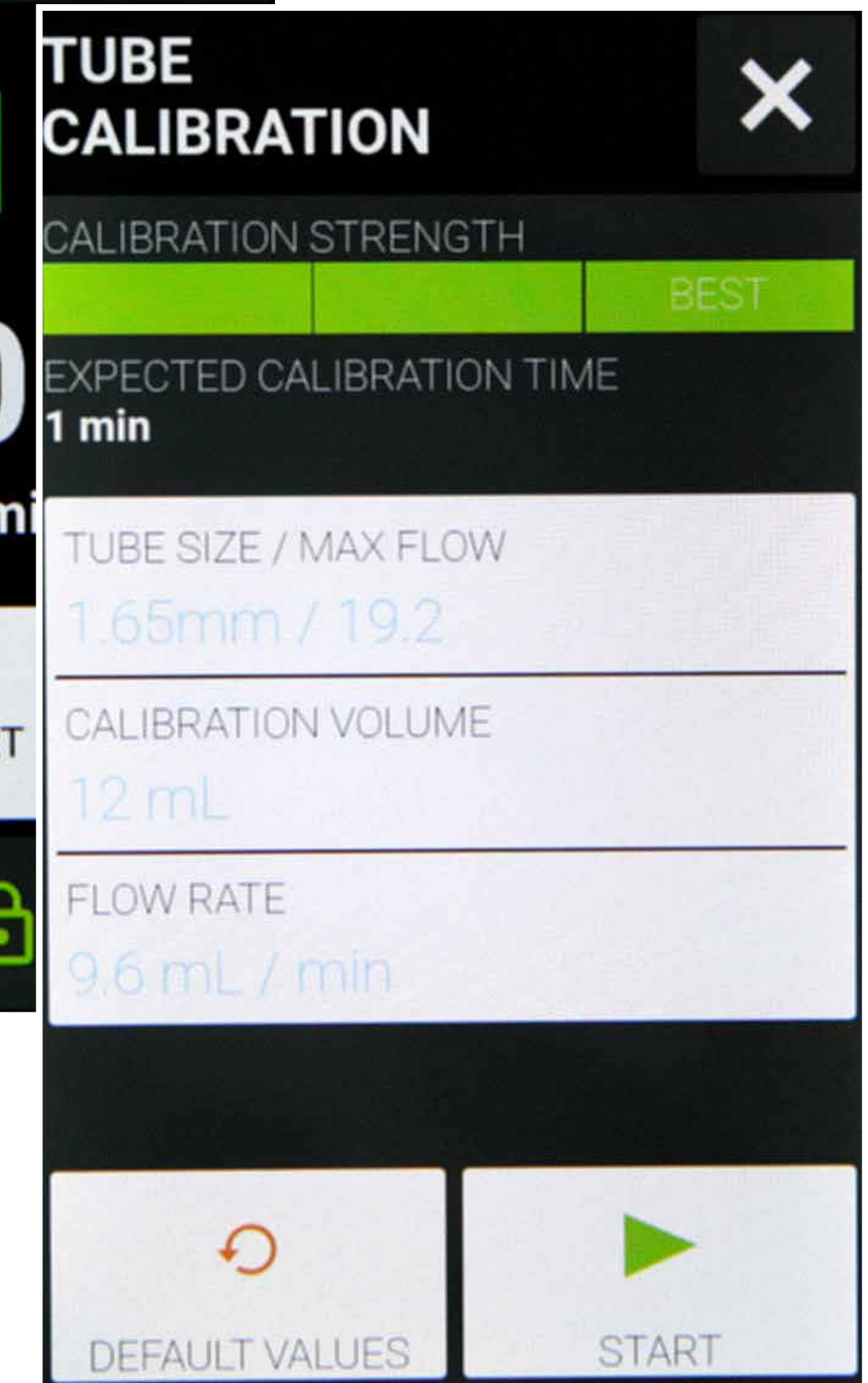

Over-The-Air Firmware Update

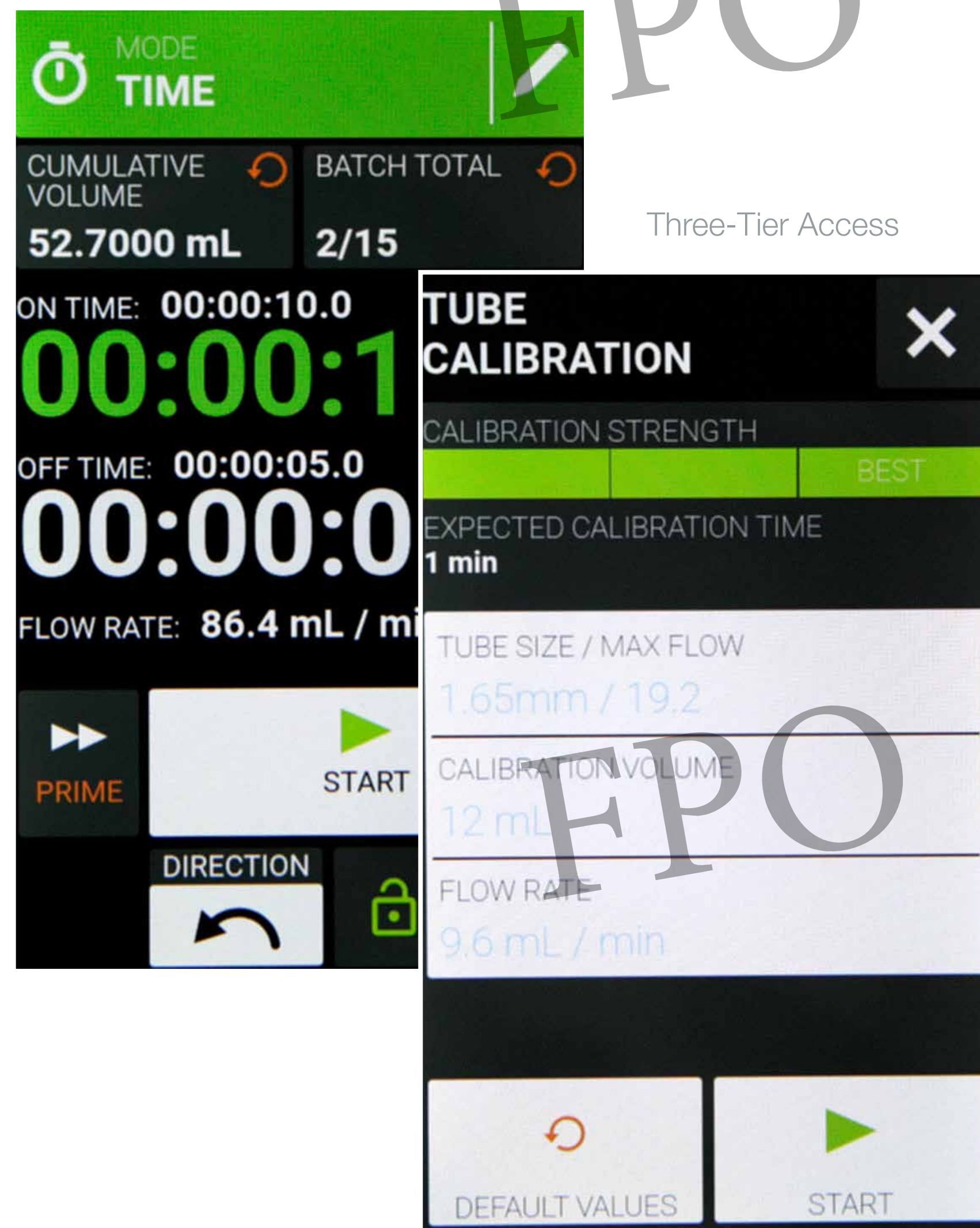

Innovative Features

- **Over-the-air firmware updates**—automatically allows you to take full advantage of future enhancements. Rest assured you will always have the latest technology on your MasterSense drive.
- **Three-tier access to pump functionality**—use this option to protect your process by controlling who can operate the pump and at what level. Offers two-factor security with user name and password required.

Connectivity Options

Get your preferred way to work with **multiple ways to operate the drive**—via touch screen, analog remote, or choose a connectivity option. When you need remote pump control and monitoring, our cloud-based MasterflexLive® technology would be the ideal choice. For a digital protocol to suit your automation needs, we offer EtherNet/IP or Profibus capability.

- **Choice of connectivity options**—MasterflexLive®, EtherNet/IP, or Profibus
- **MasterflexLive**—our secure, cloud-enabled platform for remote monitoring and control of critical processes running 24/7 without the need for onsite personnel. Push notifications alert you to changes in operating conditions and error messages. Upgrade to premium for operational data storage in our Microsoft Azure secured cloud platform that complies to both FDA 21 CFR Part 11 and EU Annex 11.

MasterflexLive is proprietary technology—patent pending.






Powerful Performance

- Broad flow range, from 0.001 to 3400 mL/min, for use **on the process floor and in the lab**
- Internet connection via Wi-Fi or standard Ethernet (RJ45) jack on rear of drive
- **Open-head sensor** technology stops the drive motor when pump head is opened or connect it to a remote leak detector to indicated ruptured filters, blown gaskets, or leaky tube connections
- Brushless, maintenance-free motor provides $\pm 0.1\%$ speed control accuracy with 6000:1 turndown ratio
- Bidirectional flow and self-priming capabilities allow for **smooth, seamless operation**


Versatility— A Pump Head for Every Application

Build the pump system for your
specific needs with ten styles of
Masterflex L/S pump head available.
**All are compatible with the
MasterSense drive.**

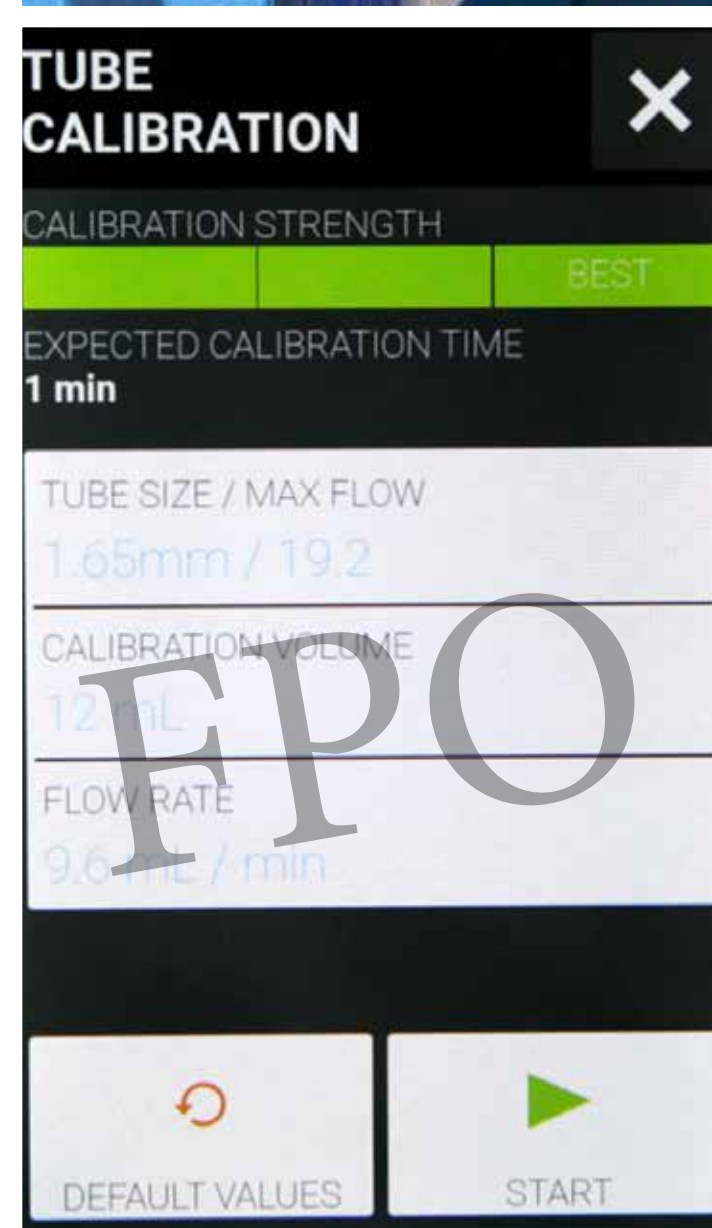


Pump Head Selection

High-Performance

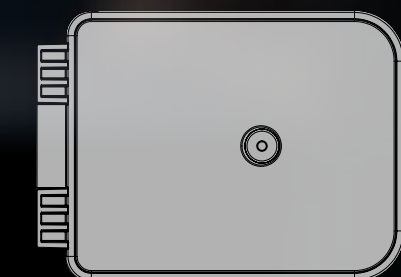


Easy-Load®

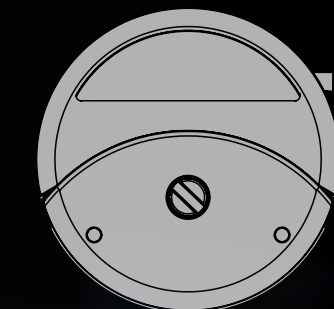


Easy-Load® 3




Easy-Load® II

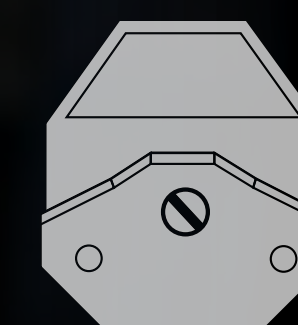


Cytoflow®

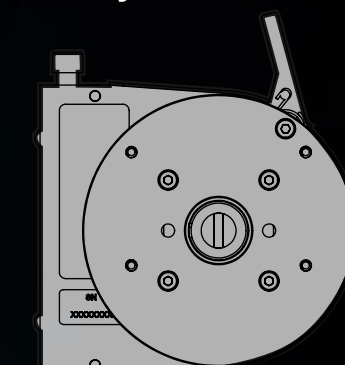


Multichannel

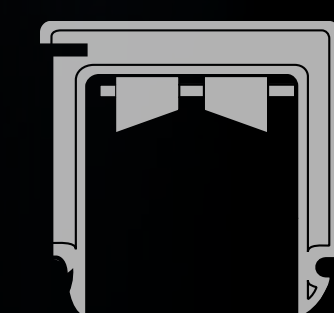


Standard

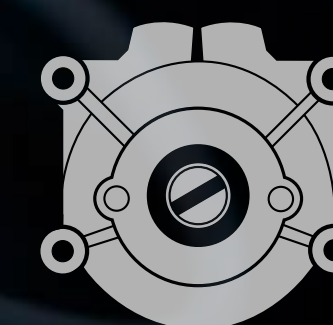

PTFE-Tubing

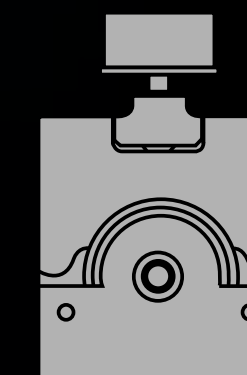

PTFE-Diaphragm

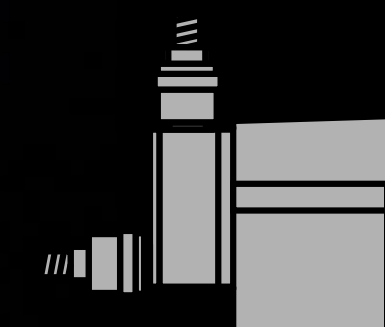

Puri-Flex™

- Heat sealable, weldable, and a great alternative to leading TPE's
- Longer pump life and lower spallation when compared to silicone or many other TPE tubings
- USP Class VI, FDA, 3A, REACH, RoHS, ADCF

Platinum-Cured Silicone

- Low protein binding, fungus resistant, phthalate and latex-free
- Accommodates a wide temperature range
- USP Class V, EP, FDA, REACH, RoHS, ADCF

Tygon® E-Food (B-44-4X)

- Smooth inner surface unaffected by commercial sanitizers
- Nontoxic; will not affect taste or odor
- FDA, EU Food, NSF, 3A, REACH, RoHS, ADCF

Viton®

- Perfect for applications where FDA and chemical compatibility is needed
- Resists corrosives, solvents, and oils at elevated temperatures
- FDA, ADCF

Masterflex size	Puri-Flex™ 25 ft (7.6 m)	Platinum-cured silicone 25 ft (7.6 m)	Tygon® E-Food (B-44-4X) 50 ft (16.2 m)	Viton® 25 ft (7.6 m)
	Item number	Item number	Item number	Item number
Precision pump tubing				
L/S 13	MK-96419-13	MK-96410-13	MK-06418-13	MK-96412-13
L/S 14	MK-96419-14	MK-96410-14	MK-06418-14	MK-96412-14
L/S 16	MK-96419-16	MK-96410-16	MK-06418-16	MK-96412-16
L/S 25	MK-96419-25	MK-96410-25	MK-06418-25	MK-96412-25
L/S 17	MK-96419-17	MK-96410-17	MK-06418-17	MK-96412-17
L/S 18	MK-96419-18	MK-96410-18	MK-06418-18	MK-96412-18
High-Performance Precision pump tubing				
L/S 15	MK-96419-15	MK-96410-15	MK-06418-15	MK-96412-15
L/S 24	MK-96419-24	MK-96410-24	MK-06418-24	MK-96412-24
L/S 35	MK-96419-35	MK-96410-35	MK-06418-35	MK-96412-35
L/S 36	MK-96419-36	MK-96410-36	MK-06418-36	MK-96412-36

Versatility—Tubing for Every Application

Optimize your demanding peristaltic pump applications when you use Masterflex pump tubing—specifically designed and manufactured for Masterflex peristaltic pumps. Select Precision tubing or the thicker-walled High-Performance Precision tubing which improves pressure generation, suction lift, ability to pump viscous fluids, and tubing life. See all 20+ tubing formulations at [Masterflex.com/tubing](https://www.masterflex.com/tubing).



Technical Specs & Ordering Information

System with Easy-Load II Pump Head



System with Easy-Load Pump Head

Flow rates in mL/min (flow rates in parentheses can be obtained only with the L/S High-Performance pump head)

rpm	Precision pump tubing						High-Performance Precision pump tubing			
	L/S 13	L/S 14	L/S 16	L/S 25	L/S 17	L/S 18	L/S 15	L/S 24	L/S 35	L/S 36
0.1 to 600	0.006 to 36	0.021 to 130	0.08 to 480	0.17 to 1000	0.28 to 1700	0.38 to 2300	0.17 to 1000 (0.18 to 1100)	0.28 to 1700 (0.30 to 1800)	0.38 to 2300 (0.43 to 2600)	0.48 to 2900 (0.58 to 3400)
0.02 to 100	0.001 to 6	0.004 to 21	0.016 to 80	0.033 to 170	0.057 to 280	0.077 to 380	0.033 to 170 (0.036 to 180)	0.057 to 280 (0.060 to 300)	0.077 to 380 (0.086 to 430)	0.096 to 480 (0.113 to 580)

Specifications

Speed control accuracy (resolution)	±0.1% (0.1 rpm at 600 rpm; 0.01 rpm at 100 rpm)
Display	5" capacitive TFT/RGB touch screen with in-plane switching (wide viewing angle); 800 x 480 pixels
External control - input	0 to 20 mA, 20 to 0 mA, 4 to 20 mA, 20 to 4 mA; 0 to 10 V, 10 to 0 V, 0 to 5 V, 5 to 0 V; Start/stop, Direction (CW/CCW), Prime via contact closure
External control - output	0 to 20 mA, 20 to 0 mA, 4 to 20 mA, 20 to 4 mA; 0 to 10 V, 10 to 0 V, 0 to 5 V, 5 to 0 V; Tach output (open collector), Motor running, Open head, Local/Remote, General alarm
Motor	BLDC with encoder; continuous duty; reversible; ¼ hp (75 W)
Open-head sensor	Yes; extension cable included
Housing materials	ABS plastic housing, painted aluminum chassis, painted aluminum bezel
Enclosure rating	IP33
Agency approvals	UL, cUL, CE
Warranty	5 years
Dimensions (L x W x H)	10¼" x 7⅝" x 8½" (26 x 19.4 x 21.6 cm)
Shipping weight	11 lb (5 kg)
Power	90 to 260 VAC, 50/60 Hz

MasterSense Drive Ordering Information

rpm	Flow range	Pumps heads accepted	With MasterflexLive	With EtherNet/IP	With Profibus
			Item number	Item number	Item number
0.1 to 600	0.006 to 3400 mL/min	2	MK-07526-16	MK-07526-10	MK-07526-12
0.02 to 100	0.001 to 580 mL/min	4	MK-07526-26	MK-07526-20	MK-07526-22

MK-07523-92 Foot switch, momentary start/stop; 6-ft (1.8-m) cable

MK-07523-94 DB25 male connector. Use to create your own cable

MK-07523-95 Cable assembly for remote control; DB25 male connector and 25-ft (7.9-m) cable with stripped wire ends

MK-07523-97 Handheld remote controller; route tubing through handle for dispensing into containers

MK-17050-01 NIST-traceable calibration with data for peristaltic pump drive

Complete Systems

Flow range	What's included			Item number
	Pump head	Tubing	Drive	
0.08 to 480 mL/min	L/S Easy-Load pump head 07514-30 with open-head sensor	Puri-Flex™ pump tubing 96419-16, L/S 16, 10 ft (3 m)	L/S MasterSense drive 07526-16, 0.1 to 600 rpm, with MasterflexLive	MK-77920-60
			L/S MasterSense drive 07526-10,0.1 to 600 rpm, with EtherNet/IP	MK-77920-40
			L/S MasterSense drive 07526-12, 0.1 to 600 rpm, with Profibus	MK-77920-50
0.28 to 1700 mL/min	L/S Easy-Load II pump head 77200-32 with open-head sensor	Puri-Flex pump tubing 96419-24, L/S 24, 10 ft (3 m)	L/S MasterSense drive 07526-16, 0.1 to 600 rpm, with MasterflexLive	MK-77920-62
			L/S MasterSense drive 07526-10,0.1 to 600 rpm, with EtherNet/IP	MK-77920-42
			L/S MasterSense drive 07526-12, 0.1 to 600 rpm, with Profibus	MK-77920-52

Masterflex is a leading, global fluidics manufacturer of proprietary products—with certifications to help you stay compliant—for your complete fluidic path. Along with our wide range of products, from pumps and tubing to flowmeters and single-use components, we are your consultant and partner to help you optimize your full fluidics workflow from research to process. And, when one of our standard offerings won't fit your application, we can help you succeed with a custom-engineered OEM pump solution to meet your exact needs and budget.

Masterflex products are known for their accuracy, durability, and ease-of-use. They can be seen in facilities around the world particularly in the pharmaceutical, biopharmaceutical, food & beverage, and chemicals industries. Masterflex products can also be used in any industries R&D and process where fluidics is required.

MASTERflex[®]

ACHROM : +Tel 32 9 380 06 47 in Belgium

Email: achrom@achrom.be /www.achrom.be

